



## ISOLATION GOWN **AAMI LEVEL 3**





Personal protective gowns are designed to create a barrier to prevent and minimize infections which may be transmitted through fluid exchange.

Resistant to liquids with its barrier feature.

Resistant to moisture and bacterial penetration.

Puncture, tear and abrasion resistant.

Free from toxic substances.

Lint and dust free.



These products meet the requirements of the European Community Regulation 2016/425 where type examination is conducted by **ASOCIACION DE INVESTIGACION DE LA INDUSTRIA TEXTIL (AITEK)** Notified Body No: **0161** and thus are **CE** marked.



**Protective Clothing Category III**

#### TYPE PB [6]-B



**EN 13034:2005+A1:2009**  
Protective clothing  
against liquid chemicals



**EN 14126:2003+AC:2004**  
Protective clothing  
against infective agents



**WHITE**



**BLUE**



### AAMI Level 3 Isolation Gown

We inform and craft the right isolation gown for the right procedure so that your choice in barrier protection becomes as intuitive as the care you dedicate to patients. AAMI Level 3 isolation gowns are intended to be used in procedures with a moderate risk of exposure.

The front closure is provided with a zipper and the zipper fly is fixed to the opposite wing with double - sided tape.

Easy to put on and take off.

Reinforcement strips on the overalls prevent leakage from the seams.



#### Velcro Tape

Velcro tape prevents opening of the collar.




#### Four Lases

Four-laced belt system insures secure fit.



#### Ribanned Manshet

Elastic rib insert at cuff maximizes fit and protection.

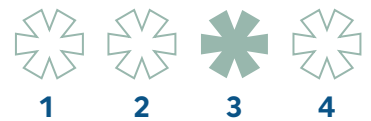

AAMI Levels determined per  
**ANSI/AAMI PB70:2012**







**SG 1389 - Type PB6-B**



0161

TYPE PB [6]-B



**Category III**


EN ISO 13688

EN 13034

EN 14126

EN 14325

EN 13274

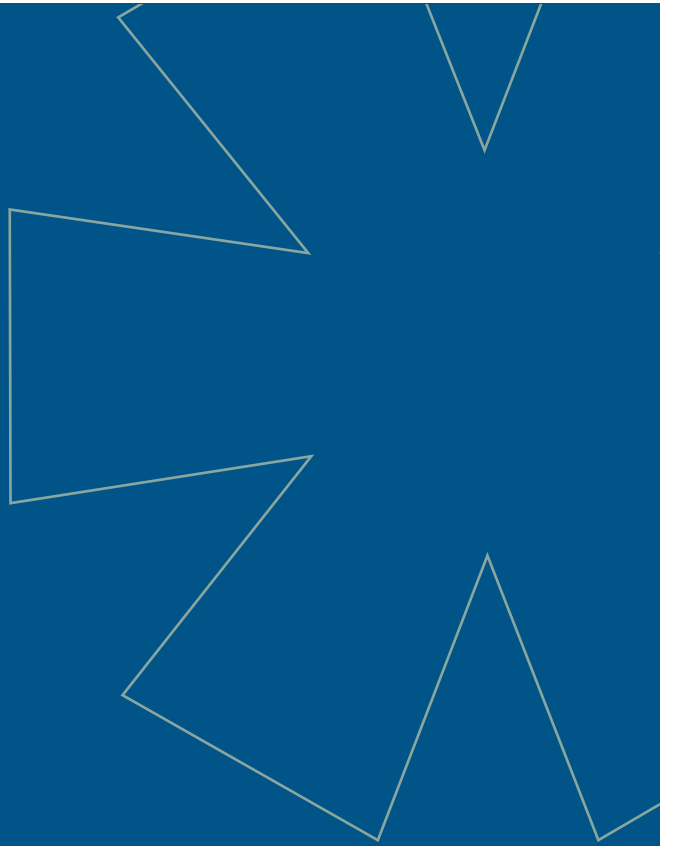

d\*lab  
BIOCLEAN

[d-lab.com](http://d-lab.com)